

Driver Software for USB Camera Module

INSTALLATION MANUAL

- Read this instruction manual carefully before using the product.
- Attempting to operate this product by methods different from those described in this instruction manual may cause serious accidents.
- Keep this instruction manual near the product for easy reference when needed.

Document Number : K46-B60051
Document Revision : 2.02 29.Jul.2014

HAMAMATSU

COPYRIGHT© HAMAMATSU PHOTONICS K.K.

Software License Terms

Please carefully review the following terms and conditions.

Hamamatsu Photonics permits the use of this software only to those companies which accepts the License Terms for using "MCD Head Products" as well as "Amprefer Board Products" ("our products"). By copying these files, or installing this software, user agrees to be bound by the License Terms. In case, the user disagrees with the complete or a part of the License Terms, he will not be able to install or use this software.

1. The purpose of this software

Considering the convenience of the users who wants to use our product easily, the software is licensed free of charge with no warranty. Please use this software at your own discretion and responsibility.

2. License

Only the user who agrees with the License Terms has the right to install and use our software application solely for purposes of implementing the measurement and control of our products.

3. Ownership of copyrights and other rights

Proprietary software and related documentation, and all other intellectual property rights will belong to us. International treaties and laws on intellectual property rights and copyright law protect this software. The user is not allowed to change or modify the attached Software or documentations.

Except for matters that are explicitly licensed under License Terms, we will not transfer or grant any rights to users. All rights on this software and documentation are reserved to us.

4. Reproduction

Users can reproduce this software for backup purposes in compliance with the terms of any provision of this license.

5. Prohibition clause

The user can not do the following things. However, in case user rents or leases or transfers our products to a third party user and the third party agrees to the License Terms of the product of our company, our company continues to permit the use of this software to a third party under this License Terms.

- 1) Sublicense, distribute or otherwise transfer the Product or any component thereof to any third party.
 - 2) Re-license or transfer the rights to use this software to the third party.
 - 3) Lease, rent or loan the Product to any unlicensed third party.
 - 4) Modify or remove part or all of this software including this License Terms and other attached annexes. The change in the file name is included in the modification.
 - 5) Reverse engineer, decompile, disassemble, or otherwise attempt to derive the source code of this Software.
-

6. Limitation of Liability

Hamamatsu Photonics will guarantee neither this software nor the attached annex including adaptability to the quality, the performance or the particular purpose at all. In any case, computer failure or damage arising from use of or inability to use this software and documentation, loss of information, all other direct and indirect damages, we will not be responsible. Moreover, not only maintenance and support but also in the event of any defect or failure, etc., we do not bear responsibility including repair and restoration of this software.

Hamamatsu Photonics can do changes without notice for improvement of this software.

Table of contents

1. Introduction	2
2. Recommended specifications	3
3. Setup	4
3.1 Setup of driver software	4
3.2 Uninstallation of driver software	7
3.3 Important note on setup of the driver software	9
3.4 Manual installation and manual update driver.	10

1

Introduction

This software installs the driver software for the USB(Universal Serial Bus) camera module control. This manual describes the installation method of the USB camera module to the PC using Windows.

Please refer to this manual while performing setup.

Recommended specifications

The following specified environment is recommended for the software.

Operating system(OS)

Microsoft® Windows 7® (x86/x64) Service Pack 1 or higher

Interface

USB2.0 compliant.

Drive device

CD-ROM drive or DVD drive.

* For install.

It explains the setup procedure of the driver software.

3.1 Setup of driver software

The installation of the driver software uses related installer software.

The required installer software is included in this CD, please install according to the following procedure.

3.1.1 Installation of driver software


The driver software is installed in the computer.

1. The installer is started.

Execute the installer file "\DriverFiles\Setup.exe" included in this CD.


2. After the installer starts, install according to the installation wizard.

The installation of the driver software begins on clicking the "Next" button on the wizard.

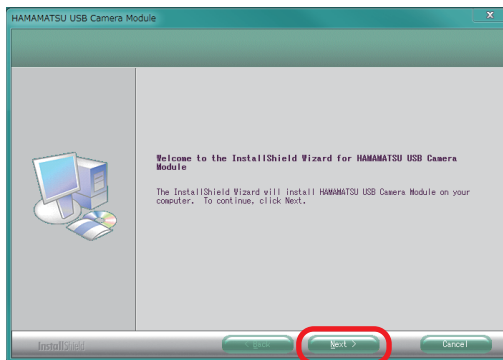


Click


3. The installation of the driver software begins on clicking the "Install" button on the wizard.

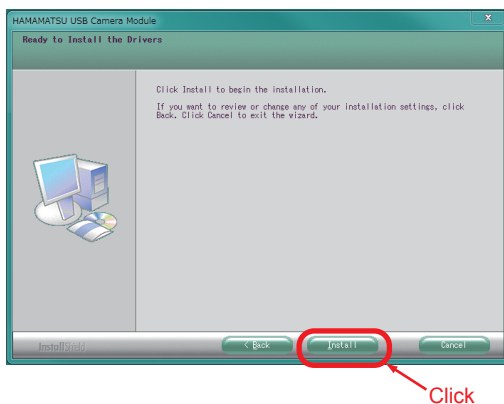

4. It takes a while for the installation procedure to complete.

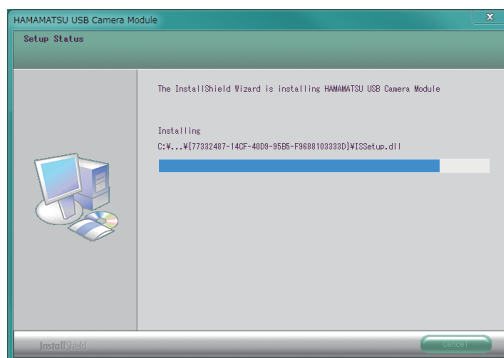


5. On successful completion, the message as shown in the figure below is displayed. Click the "Finish" Button.
The wizard is closed and successful installation of the driver software is completed.

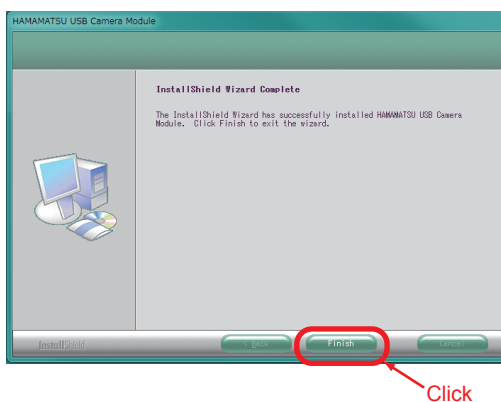


3.2 Uninstallation of driver software

The method to uninstall the driver already installed (And set up) is mentioned as follows. Uninstall the software driver using "Programs" from the "Control Panel".


1. Select "Programs" from the "Control Panel", and display the program list.
2. Select "HAMAMATSU USB Camera Module" from the displayed list.
And click the "Uninstall" button, uninstallation program start.
3. When the "OK" button on the wizard is clicked, the uninstallation of the driver software begins.

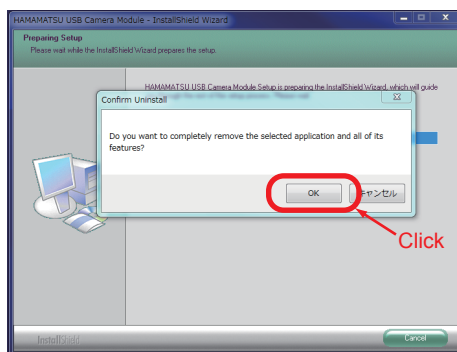

4. It takes a while until the uninstallation procedure is completed.

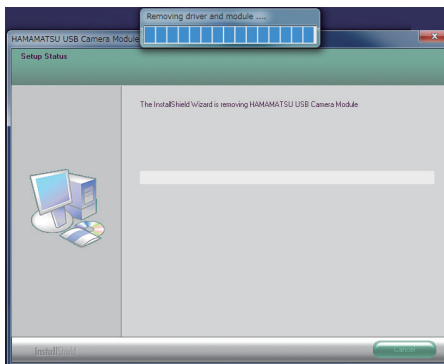

5. When the uninstallation of the driver software completes successfully, dialog is displayed, click the "Finish" button.
The wizard is closed, and uninstallation of the driver software is complete.

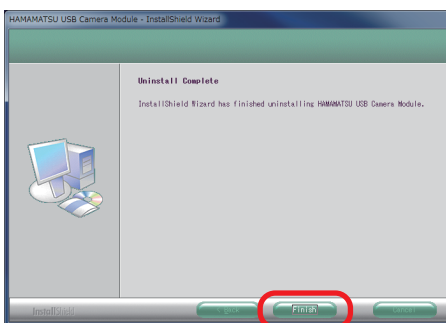

Click

3.3 Important note on setup of the driver software


When renewing the driver software to newer version or installing old version used earlier, please do so only after carrying out uninstallation of current version as per the "3.2 Uninstallation of Driver software" mentioned above.

3.4 Manual installation and manual update driver.

If you are not able to install the USB driver software, please install by manual installation. The steps of manual installation are described. In addition, if install by manual installation, the "3.2 Uninstall" is not available. (Please doing the uninstaller procedure of each OS.


3.4.1 Manual Installation (Windows7)

1. Connect USB camera module to PC with a USB cable. If the USB camera module is connected to an external power supply, connect to PC after turning on external power supply.
2. The "Welcome to the Found New Hardware Wizard" is not displayed. After connected, the product is listed on the "Device Manager", please setup the driver software by procedure of section "3.4.2 Driver Update (Windows7)".

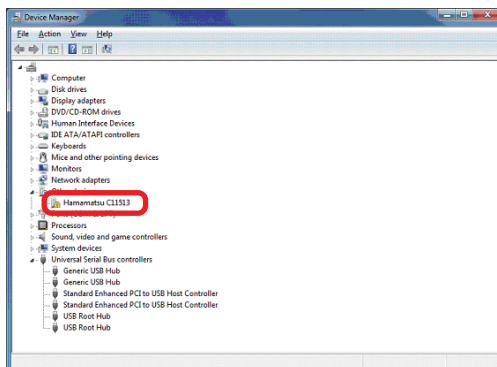


3.4.2 Driver Update (Windows7)


1. Connect USB camera module to PC with USB cable. If USB camera module is connected to an external power supply, connect to PC after turning on external power supply.
2. Insert the supplied CD-ROM into drive. Select [Start menu] → [Control panel] → [System and Security] → [System] → [Device Manager] and [Device Manager] window will be displayed. After displaying [Device Manager] window, USB camera module will be displayed by expanding [Universal Serial Bus controllers].

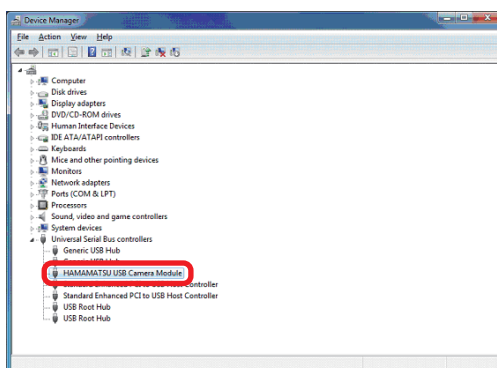


3. Right Click [HAMAMATSU USB Camera Module] and select [Update Driver Software].

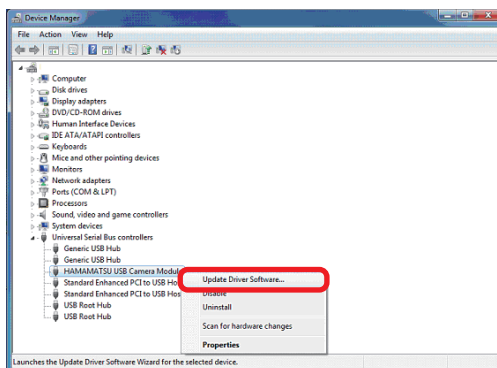


4. When [Update Driver Software] is selected, [Update Driver Software] window will be displayed. Select [Browse my computer...].

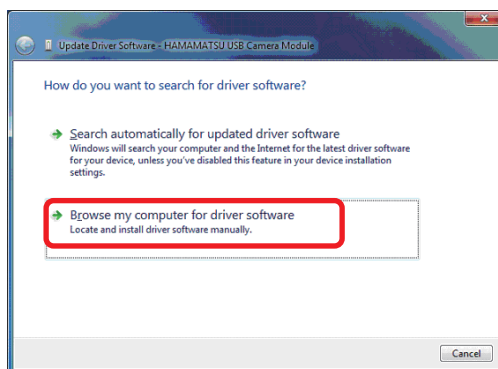


5. In this screen, the driver file contained in the CD-ROM that was set into the drive will be searched. Please click [Browse] button and from the dialog box that will be displayed, select the appropriate OS folder from the [DriverFiles] folder on the CD-ROM. After selecting the folder, click [OK] button and then click [Next] button on this screen.

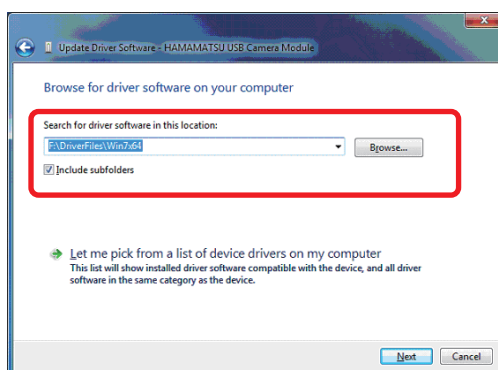

6. After installation is started, the following screen will be displayed.

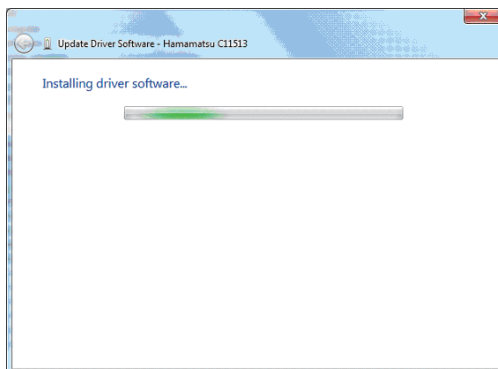

7. When installation is finished, the next screen will be displayed.
Click [Close] button.

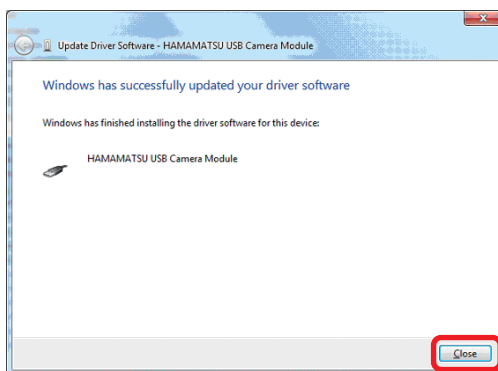

Click

8. After driver is installed, select [Start Menu] → [Control Panel] → [System and Security] → [System] → [Device Manager] and [Device Manager] window will be displayed. After displaying [Device Manager] window, please confirm that the USB camera module is displayed by expanding [Universal Serial Bus controllers].

Document History

Date	Document Revision	Contents
07.Aug.2009	1.00	First Printing
08.Oct.2009	1.01	Corrected misspelling.
28.Jun.2012	2.00	<ul style="list-style-type: none">• Delete the description of Wndows2000, and added the description of Windows7.• Added the description of manual setup.
20.Sep.2013	2.01	Add list of offices on the back cover.
29.Jul.2014	2.02	Edited by the end of support of WindowsXP.

USB Camera Module Control Driver Software Installation Manual

Manufacturer

Japan: Hamamatsu Photonics K.K., Solid State Division

1126-1 Ichino-cho, Higashi-ku, Hamamatsu City, 435-8558 Japan, Telephone: (81)53-434-3311, Fax: (81)53-434-5184
E-mail: export@sys.hpk.co.jp

U.S.A. and Canada: Hamamatsu Corporation, Systems Division

360 Foothill Road, Bridgewater, N.J. 08807-0910, U.S.A. Telephone: (1)908-231-1116, Fax: (1)908-231-0852
E-mail: usa@hamamatsu.com

Germany: Hamamatsu Photonics Deutschland GmbH

Arzbergerstr. 10, D-82211 Herrsching am Ammersee, Germany Telephone: (49)8152-375-0, Fax: (49)8152-265-8
E-mail: info@hamamatsu.de

France: Hamamatsu Photonics France S.A.R.L.

19, Rue du Saule Trapu, Parc du Moulin de Massy, 91882 Massy Cedex, France Telephone: (33)1 69 53 71 00, Fax: (33)1 69 53 71 10
E-mail: infos@hamamatsu.fr

United Kingdom: Hamamatsu Photonics UK Limited

2 Howard Court, 10 Tewin Road, Welwyn Garden City Hertfordshire AL7 1BW, United Kingdom Telephone: (44)1707-294888, Fax: (44)1707-325777 E-mail: info@hamamatsu.co.uk

North Europe: Hamamatsu Photonics Norden AB

Smidesvagan 12, SE-171 41 Solna, Sweden Telephone: (46)8-509-031-00, Fax: (46)8-509-031-01
E-mail: info@hamamatsu.se

Italy: Hamamatsu Photonics Italia S.R.L.

Strada della Moia, 1/E 20020 Arese (Milano), Italy Telephone: (39)02-935 81 733, Fax: (39)02-935 81 741

WEB SITE : <http://www.hamamatsu.com/>

Document Number : K46-B60051

Document Revision 2.02 29.Jul.2014

Information in this document is subject to change without notice.

Microsoft®, Windows®, Windows XP®, Windows 7® are Microsoft Corporation trademarks or registered trademarks. Intel®, Pentium® are Intel Corporation trademarks or registered trademarks.